



GenWay Biotech, Inc.

6777 Nancy Ridge Drive

San Diego, CA 92121

Phone: 858.458.0866

Fax: 858.458.0833

Email: techline@genwaybio.com

<http://www.genwaybio.com>

Human BDNF ELISA Kit

Catalog No. **GWB-ZZD022**

Size 96T(8×12 divisible strips)

For quantitative detection of human BDNF in cell culture supernates, serum and plasma(heparin, EDTA, citrate).


Typical Data Obtained from Human BDNF

(TMB reaction incubate at 37°C for 25-30 min)

Concentration(pg/ml)	0.0	31.2	62.5	125	250	500	1000	2000
O.D	0.034	0.062	0.104	0.215	0.462	1.003	1.840	2.565

Typical Human BDNF ELISA Kit Standard Curve

This standard curve was generated at GenWay Biotech, Inc. for demonstration purpose only. A standard curve must be run with each assay.

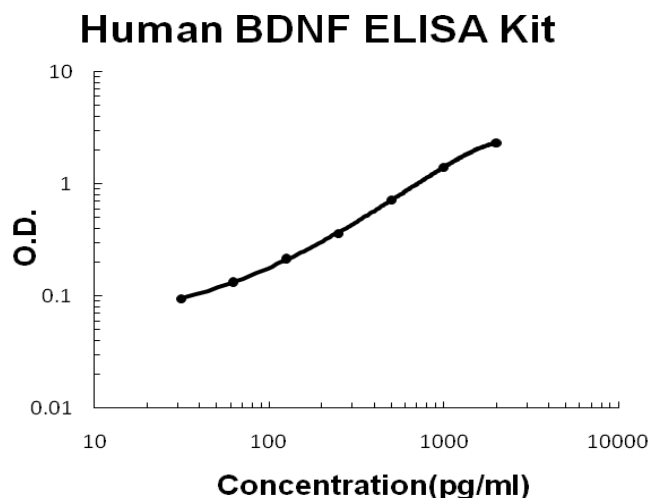

Range	31.2pg/ml-2000pg/ml
Sensitivity	< 2pg/ml
Specificity	Natural and recombinant human BDNF
Cross-reactivity	No detectable cross-reactivity with other relevant proteins

Storage

Store at 4°C for 6 months, at -20°C for 12 months. Avoid multiple freeze-thaw cycles (Shipped with wet ice.)

FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC AND CLINICAL USE.

Product Information Sheet

Precision

Intra-Assay Precision (Precision within an assay) Three samples of known concentration were tested on one plate to assess intra-assay precision.

Inter-Assay Precision (Precision between assays) Three samples of known concentration were tested in separate assays to assess inter-assay precision.

Sample	Intra-Assay Precision			Inter-Assay Precision		
	1	2	3	1	2	3
n	16	16	16	24	24	24
Mean(pg/ml)	113	425	1278	197	594	1306
Standard deviation	3.96	17.85	62.62	14.18	44.55	103.17
CV(%)	3.5	4.2	4.9	7.2	7.5	7.9

Principle

GenWay Biotech, Inc.'s human BDNF ELISA Kit was based on standard sandwich enzyme-linked immune-sorbent assay technology. A monoclonal antibody from mouse specific for BDNF has been precoated onto 96-well plates. Standards(sf21, H129-R247) and test samples are added to the wells, a biotinylated detection polyclonal antibody from goat specific for BDNF is added subsequently and then followed by washing with PBS or TBS buffer. Avidin-Biotin-Peroxidase Complex was added and unbound conjugates were washed away with PBS or TBS buffer. HRP substrate TMB was used to visualize HRP enzymatic reaction. TMB was catalyzed by HRP to produce a blue color product that changed into yellow after adding acidic stop solution. The density of yellow is proportional to the human BDNF amount of sample captured in plate.

Kit Components

Description	Quantity
96-well plate precoated with anti- human BDNF antibody	1
Lyophilized recombinant human BDNF standard	10ng/tubex2
Biotinylated anti- human BDNF antibody	130µl(dilution 1:100)
Avidin-Biotin-Peroxidase Complex(ABC)	130µl(dilution 1:100)
Sample diluent buffer	30 ml
Antibody diluent buffer	12ml
ABC diluent buffer	12ml
TMB color developing agent	10ml
TMB stop solution	10ml

Material Required But Not Provided

1. Microplate reader in standard size.
2. Automated plate washer.
3. Adjustable pipettes and pipette tips. Multichannel pipettes are recommended in the condition of large amount of samples in the detection.
4. Clean tubes and Eppendorf tubes.
5. Washing buffer (neutral PBS or TBS).
 - Preparation of 0.01M **TBS**: Add 1.2g Tris, 8.5g NaCl; 450µl of purified acetic acid or 700µl of concentrated hydrochloric acid to 1000ml H₂O and adjust pH to 7.2-7.6. Finally, adjust the total volume to 1L.
 - Preparation of 0.01 M **PBS**: Add 8.5g sodium chloride, 1.4g Na₂ HPO₄ and 0.2g NaH₂ PO₄ to 1000ml distilled water and adjust pH to 7.2-7.6. Finally, adjust the total volume to 1L.

Notice for Application of Kit

1. To inspect the validity of experiment operation and the appropriateness of sample dilution proportion, pilot experiment using standards and a small number of samples is recommended.
2. The TMB Color Developing agent is colorless and transparent before using, contact us freely if it is not the case.
3. Before using the Kit, spin tubes and bring down all components to the bottom of tubes.
4. Duplicate well assay is recommended for both standard and sample testing.
5. Don't let 96-well plate dry, for dry plate will inactivate active components on plate.
6. Don't reuse tips and tubes to avoid crosscontamination.
7. Avoid using the reagents from different batches together.
8. In order to avoid marginal effect of plate incubation due to temperature difference (reaction may be stronger in the marginal wells), it is suggested that the diluted ABC and TMB solution will be pre-warmed in 37°C for 30 min before using.

Preparation

1. Sample Preparation and Storage

Store samples to be assayed within 24 hours at 2-8°C. For long-term storage, aliquot and freeze samples at -20°C.

Avoid repeated freeze-thaw cycles.

- **Cell culture supernates:** Remove particulates by centrifugation, assay immediately or aliquot and store samples at -20°C.
- **Serum:** Allow the serum to clot in a serum separator tube (about 4 hours) at room temperature. Centrifuge at approximately 1000 X g for 15 min. Analyze the serum immediately or aliquot and store samples at -20°C.
- **Plasma:** Collect plasma using **heparin, EDTA** or **citrate** as an anticoagulant. Centrifuge for 15 min at 1500 x g within 30 min of collection. Assay immediately or aliquot and store samples at -20°C.

2. Sample Dilution Guideline

The user needs to estimate the concentration of the target protein in the sample and select a proper dilution factor so that the diluted target protein concentration falls near the middle of the linear regime in the standard curve. Dilute the sample using the provided diluent buffer. The following is a guideline for sample dilution. Several trials may be necessary in practice. **The sample must be well mixed with the diluents buffer.**

- **High target protein concentration (20-200ng/ml).** The working dilution is 1:100. i.e. Add 1µl sample into 99 µl sample diluent buffer.
- **Medium target protein concentration (2-20ng/ml).** The working dilution is 1:10. i.e. Add 10µl sample into 90 µl sample diluent buffer.
- **Low target protein concentration (31.2-2000pg/ml).** The working dilution is 1:2. i.e. Add 50µl sample to 50 µl sample diluent buffer.
- **Very Low target protein concentration \leq (31.2pg/ml).** No dilution necessary, or the working dilution is 1:2.

3. Reagent Preparation and Storage

A. Reconstitution of the human BDNF standard: BDNF standard solution should be prepared no more than 2 hours prior to the experiment. Two tubes of BDNF standard (10ng per tube) are included in each kit. Use one tube for each experiment.

- a. 10,000pg/ml of human BDNF standard solution: Add 1ml sample diluent buffer into one tube, keep the tube at room temperature for 10 min and mix thoroughly.
- b. 2000pg/ml of human BDNF standard solution: Add 0.2ml of the above 10ng/ml BDNF standard solution into 0.8 ml sample diluent buffer and mix thoroughly.

FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC AND CLINICAL USE.

Product Information Sheet

- c. 1000pg/ml→31.2pg/ml of human BDNF standard solutions: Label 6 Eppendorf tubes with 1000pg/ml, 500pg/ml, 250pg/ml, 125pg/ml, 62.5pg/ml, 31.2pg/ml respectively. Aliquot 0.3ml of the sample diluent buffer into each tube. Add 0.3ml of the above 2000pg/ml BDNF standard solution into 1st tube and mix. Transfer 0.3ml from 1st tube to 2nd tube and mix. Transfer 0.3ml from 2nd tube to 3rd tube and mix, and so on.

Note: The standard solutions are best used within 2 hours. The 10ng/ml standard solution should be stored at 4°C for up to 12 hours, or at -20°C for up to 48 hours. Avoid repeated freeze-thaw cycles.

- B. Preparation of biotinylated anti-human BDNF antibody working solution: The solution should be prepared no more than 2 hours prior to the experiment.
 - a. The total volume should be: 0.1ml/well x (the number of wells). (Allowing 0.1-0.2 ml more than total volume)
 - b. Biotinylated anti-human BDNF antibody should be diluted in 1:100 with the antibody diluent buffer and mixed thoroughly. (i.e. Add 1µl Biotinylated anti-human BDNF antibody to 99µl antibody diluent buffer.)
- C. Preparation of Avidin-Biotin-Peroxidase Complex (ABC) working solution: The solution should be prepared no more than 1 hour prior to the experiment.
 - a. The total volume should be: 0.1ml/well x (the number of wells). (Allowing 0.1-0.2 ml more than total volume)
 - b. Avidin- Biotin-Peroxidase Complex (ABC) should be diluted in 1:100 with the ABC dilution buffer and mixed thoroughly. (i.e. Add 1µl ABC to 99µl ABC diluent buffer.)

Assay Procedure

The ABC working solution and TMB color developing agent must be kept warm at 37°C for 30 min before use. When diluting samples and reagents, they must be mixed completely and evenly. Standard BDNF detection curve should be prepared for each experiment. The user will decide sample dilution fold by crude estimation of BDNF amount in samples.

1. Aliquot 0.1ml per well of the 2000pg/ml, 1000pg/ml, 500pg/ml, 250pg/ml, 125pg/ml, 62.5pg/ml, 31.2pg/ml human BDNF standard solutions into the precoated 96-well plate. Add 0.1ml of the sample diluent buffer into the control well (Zero well). Add 0.1ml of each properly diluted sample of human cell culture supernates, serum or plasma (heparin, EDTA, citrate) to each empty well. **See “Sample Dilution Guideline” above for details.** It is recommended that each human BDNF standard solution and each sample be measured in duplicate.
2. Seal the plate with the cover and incubate at 37°C for 90 min.
3. Remove the cover, discard plate content, and blot the plate onto paper towels or other absorbent material. Do NOT let the wells completely dry at anytime.
4. Add 0.1ml of biotinylated anti-human BDNF antibody working solution into each well and incubate the plate at 37°C for 60 min.
5. Wash plate 3 times with 0.01M TBS or 0.01M PBS, and each time let washing buffer stay in the wells for 1 min. Discard the washing buffer and blot the plate onto paper towels or other absorbent material. (**Plate Washing Method:** Discard the solution in the plate without touching the side walls. Blot the plate onto paper towels or other absorbent material. Soak each well with at least 0.3 ml PBS or TBS buffer for 1~2 minutes. Repeat this process two additional times for a total of THREE washes. Note: For automated washing, aspirate all wells and wash THREE times with PBS or TBS buffer, overfilling wells with PBS or TBS buffer. Blot the plate onto paper towels or other absorbent material.)
6. Add 0.1ml of prepared ABC working solution into each well and incubate the plate at 37°C for 30 min.
7. Wash plate 5 times with 0.01M TBS or 0.01M PBS, and each time let washing buffer stay in the wells for 1-2 min. Discard the washing buffer and blot the plate onto paper towels or other absorbent material.

(See Step 5 for plate washing method).

8. Add 90µl of prepared TMB color developing agent into each well and incubate plate at 37°C in dark for 25-30 min (**Note:** For reference only, the optimal incubation time should be determined by end user. And the shades of blue can be seen in the wells with the four most concentrated human BDNF standard solutions; the other wells show no obvious color).
9. Add 0.1ml of prepared TMB stop solution into each well. The color changes into yellow immediately.
10. Read the O.D. absorbance at 450nm in a microplate reader within 30 min after adding the stop solution.

For calculation, (the relative O.D.₄₅₀) = (the O.D.₄₅₀ of each well) – (the O.D.₄₅₀ of Zero well). The standard curve can be plotted as the relative O.D.₄₅₀ of each standard solution (Y) vs. the respective concentration of the standard solution (X). The human BDNF concentration of the samples can be interpolated from the standard curve. **Note:** if the samples measured were diluted, multiply the dilution factor to the concentrations from interpolation to obtain the concentration before dilution.

Summary

1. Add samples and standards and incubate the plate at 37°C for 90 min. Do not wash.
2. Add biotinylated antibodies and incubate the plate at 37°C for 60 min. Wash plate 3 times with 0.01M TBS.
3. Add ABC working solution and incubate the plate at 37°C for 30 min. Wash plate 5 times with 0.01M TBS.
4. Add TMB color developing agent and incubate the plate at 37°C in dark for 25-30 min.
5. Add TMB stop solution and read.

Background

Brain-derived neurotrophic factor (BDNF) is a prosurvival factor induced by cortical neurons that is necessary for survival of striatal neurons in the brain. It is a secreted protein with the molecular weight of 27.8kDa, consisting of 247 amino acids. It is known to promote neuronal survival and differentiation. BDNF shares substantial amino acid sequence identity with nerve growth factor (NGF). BDNF and neurotrophin-3 (NT-3) are two recently cloned neurotrophic factors that are homologous to NGF. mRNA products of the BDNF and NT-3 genes are detected in the adult human brain, suggesting that these proteins are involved in the maintenance of the adult nervous system.¹ BDNF and other neurotrophins are critically involved in long-term potentiation (LTP). BDNF-mediated LTP is induced postsynaptically.² BDNF has trophic effects on serotonergic (5-HT) neurons in the central nervous system.³ BDNF has an essential maintenance function in the regulation of anxiety-related behavior and in food intake through central mediators in both the basal and fasted state.⁴ It plays a role in treating breathing disorders such as respiratory insufficiency after spinal injury.⁵ The mature form of BDNF is identical in all mammals examined, and the gene encoding human BDNF to chromosome 11, bandp13.⁶

Reference

1. Jones, K. R.; Reichardt, L. F. Molecular cloning of a human gene that is a member of the nerve growth factor family. *Proc. Nat. Acad. Sci.* 87: 8060-8064, 1990.
2. Kovalchuk, Y.; Hanse, E.; Kafitz, K. W.; Konnerth, A. Postsynaptic induction of BDNF-mediated long-term potentiation. *Science* 295: 1729-1734, 2002.
3. Lyons, W. E.; Mamounas, L. A.; Ricaurte, G. A; Coppola, V.; Reid, S. W.; Bora, S. H.; Wihler, C.; Koliatsos, V. E.; Tessarollo, L. Brain-derived neurotrophic factor-deficient mice develop aggressiveness and hyperphagia in conjunction with brain serotonergic abnormalities. *Proc. Nat. Acad. Sci.* 96: 15239-15244, 1999.
4. Rios, M.; Fan, G.; Fekete, C.; Kelly, J.; Bates, B.; Kuehn, R.; Lechan, R. M.; Jaenisch, R. Conditional deletion of brain-derived neurotrophic factor in the postnatal brain leads to obesity and hyperactivity. *Molec. Endocr.* 15: 1748-1757, 2001.
5. Baker-Herman, T. L.; Fuller, D. D.; Bavis, R. W.; Zabka, A. G.; Golder, F. J.; Doperalski, N. J.; Johnson, R. A.; Watters, J. J.; Mitchell, G. S. BDNF is necessary and sufficient for spinal respiratory plasticity following intermittent hypoxia. *Nature Neurosci.* 7: 48-55, 2004.
6. Maisonpierre, P. C.; Le Beau, M. M.; Espinosa, R., III; Ip, N. Y.; Belluscio, L.; de la Monte, S. M.; Squinto, S.; Furth, M. E.; Yancopoulos, G. D. Human and rat brain-derived neurotrophic factor and neurotrophin-3: gene structures, distributions and chromosomal localizations. *Genomics* 10: 558-568, 1991.